

Information COD Cuvette Test

**New:
ST-COD
according
to ISO**

COD analysis now Approved, practical, affordable: ST-COD – conforming to ISO 15705

The advantages


- **Approved method:** The HACH COD cuvette test is based on the same principles as the reference method
- **Comparable results** to the reference method
- **Simple, time-saving operation:** All reagents are ready to use, in the right amounts, in the cuvette
- **High level of operational safety** as the test is a closed system
- **Environmental protection** thanks to minimal use of chemicals
- **Quality** made in Germany

Compact COD


The sum parameter chemical oxygen demand (COD) is a measure of the organic pollution of wastewater. It is the most frequently mentioned parameter in European legislation. It serves as a planning assessment basis and is used to monitor the degradation efficiency of wastewater plants.

With the Small Tube COD (ST-COD) HACH offers the right COD solution with the exact ISO 15705 formulation.


Test procedure

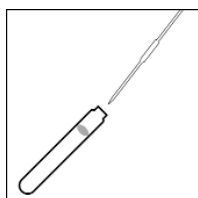


Pipette sample into cuvette

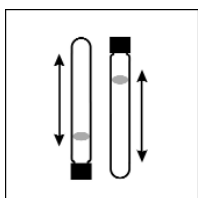

Invert cuvettes carefully

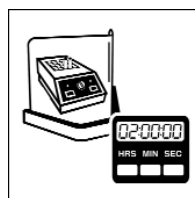

Heat cuvettes for 2 h at 150 °C

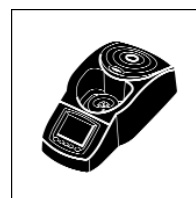

Evaluate cuvettes with photometer

Measuring ranges + applications

For industrial water, e.g. in the food, beverages and paper sectors

→ **200–15,000 mg/l**

For industrial wastewater and inflows of municipal sewage treatment plants

→ **0–1,000 mg/l (ST-COD)**

→ **20–1,500 mg/l**

For outflows of industrial and municipal sewage treatment plants

→ **0–150 mg/l (ST-COD)**

→ **3–150 mg/l**

For sewage treatment plant outflows, surface water and other water with low organic pollution:

→ **1–40 mg/l**

Information

COD Cuvette Test


The facts about the COD test speak louder than words!



Acceptance

The decisive criterion for the use of alternative test methods is their comparability with the reference method. The graph below shows a comparison of the regular HACH cuvette test and the reference - method based on analyses of 20 wastewater samples of municipal and industrial origin.

Conclusion: The correlation between the results of both methods is very good over the whole measuring range.

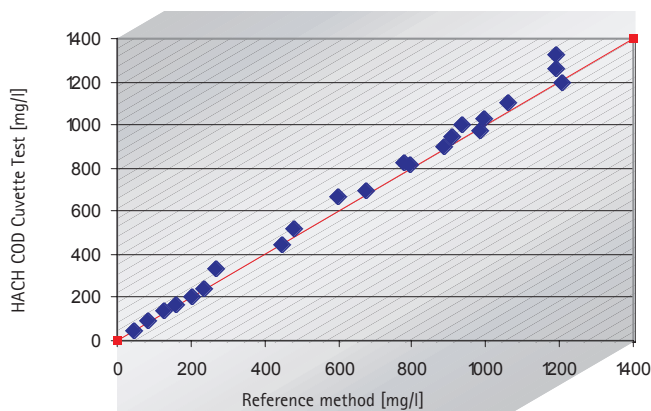

Environmental protection + safety

From an ecological point of view the COD cuvette test is far superior to the reference method. The amount of environmentally harmful chemicals it uses - especially potassium dichromate, sulphuric acid, silver and mercury - is 90 % less than the reference cuvette (see the illustration above). The risk to the user through direct contact with the chemicals is largely excluded by the closed analysis system.

Order information

Art. no.	Measuring Range
24159-51	200-15,000 mg/l
21259-51	20-1,500 mg/l
21258-51	3-150 mg/l

24158-51* 1-40 mg/l

All tests with 25 analysis cuvettes/package

* Evaluation only with DR/2500 and DR/4000

ST-COD according to ISO 15705

HCT 190** 0-1,000 mg/l

HCT 191** 0-150 mg/l

All tests with 24 analysis cuvettes/package

** Evaluation with DR/2010, DR/2500, DR/2400 and DR/4000

DR. BRUNO LANGE
GMBH & CO. KG
Willstätterstraße 11
D-40549 Düsseldorf
Tel. +49(0)211 5288-0
Fax +49(0)211 5288-143
kundenservice@drlange.de
www.drlange.de

HACH LANGE LTD
Lennox Road
Basingstoke
Hampshire, RG22 4AP
Tel. +44 (0)1256 333 403
Fax +44 (0)1256 330 724
info@drlange.co.uk
www.hach-lange.co.uk

Hyxo Oy

P.O. Box 16 (Palokorvenkatu 2)
FI-04261 Kerava, Finland
Tel. +358 10 417 4500
Fax +358 10 417 4501
hyxo@hyxo.fi • www.hyxo.com

