



SS-245

SS-325

WITH SAFETY INTERLOCK SYSTEM



Features:-

Safety Features:

- Lid Interlock System

The lid interlock system works when the electrical power is on. The lid interlock system is locked when the chamber pressure is over 0kPa.

- Stainless Steel Construction

The chamber and lid are made of stainless steel for higher mechanical strength and resistance to corrosion. The heater and baskets are also made of stainless steel.

- Ground Leakage Breaker

The power switch also acts as a circuit breaker. By breaking the power line upon detection of an excessive leakage to the ground, the possibility of electrical shock is eliminated.

- Pressure Release Valve

This valve opens to reduce chamber pressure in case of over-pressure.

Top-Loading Design:

The Top Loading TOMY Autoclaves optimize use of chamber space. At conventional horizontal-loading autoclaves, rectangular baskets are placed into a cylindrical chamber. Those situation loses 36% of the chamber volume. With a Top-Loading design, the Autoclave is installed on the floor, eliminating the need for a stand, a shelf, or valuable benchtop space for installation.

SS-245 and SS-325 do not require water and steam connection for flexible installation at anywhere, SS-325 has casters and SS-245 has rubber feet for more portability.



Microprocessor Control:

The sterilization temperature and time are controlled by a microprocessor for accurate control and reproducibility. Once sterilization temperature and time have set, data remains in memory even after power has been shut. Temperature and time displays are very clear and large for easy read-out and time display displays error warning when it happens also.

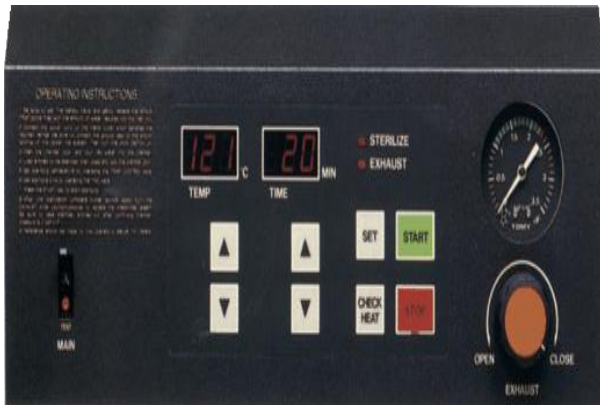


Image on Left: Operating Panel.
Image on Right: Error Indicator Display



Self-Contained Waste Bottle:

Steam vents into a waste bottle located in front of the autoclave, thus eliminating the need for a drainage hose connection and the problems associated with it. The water level of the bottle can be easily confirmed so that the water can be discarded before it overflows. Venting the chamber steam can be completed quickly or slowly over a period of 15 to 20 minutes.



SS-245

SS-325

Chamber Capacity	22 liters	53 liters
Operating Temperature	60°C to 123°C	60°C to 132°C
Operating Pressure	0 to 1.2 kgf/cm ² (0 to 17.1 psi)	0 to 1.9 kgf/cm ² (0 to 27.1 psi)
Temperature Control	Digital, Microprocessor Control	Digital, Microprocessor Control
Temperature Display Range	Digital, 0 - 180 °C	Digital, 0 - 180 °C
Pressure Display Range	0 to 3.5 kgf/cm ² (0 to 50 psi)	
Maximum Operating Pressure	1.5 kgf/cm ² (21.3 psi)	2.2 kgf/cm ² (31.2 psi)
Sterilization Heat Source	Stainless steel electric heater	Stainless steel electric heater
Safety Device	1. Lid Interlock System 2. Over Temperature System 3. Chamber Overheat Sensor 4. Overpressure Detection System 5. Insufficient Water Detection System 6. Ground LEakage Breakder Detection System	
Width	410 mm (15-3/4 in.)	500 mm (19-5/16 in.)
Depth	460 mm (18-1/8 in.)	560 mm (22-1/16 in.)
Height	905 mm (35-1/4 in.)	1090 mm (42-7/5 in.)
Chamber Diameter	248 mm (9-3/4 in.)	325 mm (12-3/4 in.)
Chamber Height	450 mm (17-11/16 in.)	635 mm (25 in.)
Power (Other voltages avail. upon request)	220 V, 50/60 Hz, 8 amps	220 V, 50/60 Hz, 10 amps